

Purair[®] FLEX

Portable, Cost Effective Flexible Film Isolator

"Innovative Design for Forensic and Biological Investigations."



— Model FLEX-24

Flexible Containment For:
Forensic Evidence Investigation;
Biological Contamination Containment and
Repair or Manipulation of Contaminated Equipment



THE AIR SCIENCE PERFORMANCE ADVANTAGE

Air Science provides containment solutions to meet any analytical need. The Purair Flex provides numerous performance advantages over the competition, including:

More Volume.
The Purair Flex Film Isolator features an innovative curved film design that provides more working volume than any other Glove Bag on the market.

Versatile Tapered Sleeve / Glove Options.
Unlike other glove bags, the double O-ring design on the standard polyurethane cuffs allow users to quickly and easily change gloves to meet a variety of dexterity needs.

Easy Setup.
Semi-rigid support rods make the Purair Flex easy to setup up and provide increased stability, even if the bag is not inflated to full pressure.


Safety Options.
HEPA filter availability, Bag-in / Bag-out Port, and optional Nitrogen purge inlet connections offer the reliability and safety today's analysts need.




PRODUCT FEATURES:

- A. Tag Sleeve on Rear of Unit:** The tag sleeve on the rear of the unit serves as a convenient method to label dangerous or sensitive contents to keep users or materials safe.
- B. Suction Cup Feet on Bottom of Base:** The Purair Flex comes equipped with suction cup feet to secure the flexible polypropylene work area and increase stability of the unit.
- C. Angled Zipper Entry.**
- D. Polyurethane Tapered Sleeves with Double O-ring Cuffs.**


FLEX-24

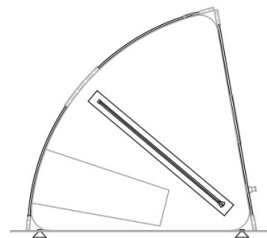

Side View

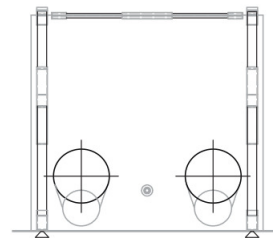

FLEX-24

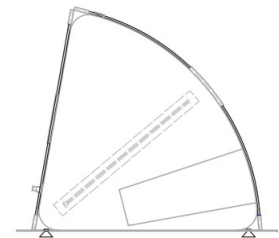

Side View

MODEL	DIMENSIONS	WEIGHT (lbs)
	External (W x D x H)	Net
Purair Flex™		
FLEX-24	23.6" x 24.6" x 23.6" 600 x 625 x 600 mm	< 10 lbs

PRODUCT SPECIFICATIONS

Purair Flex™ Model	FLEX-24
Construction Bag	Polyurethane film, Clear
Construction Base	PVC film, Blue
Sleeves	Tapered Anti-Static Polyurethane film, 450 mm long, Frosted
Gloves	Standard PVC Size XL, mounted on Double O-Ring Cuffs, replaceable
Support Frame	PVC Strips, white
Zipper	YKK brand, air pressure resistant
Pressure Test	330pa
Breather HEPA Filter, side wall mounted	Optional
Nitrogen Gas Inlet Port	Optional
Bag-in / Bag-out port (200mm diameter) with 6M liner	Optional
Tag Sleeve, rear mounted	Standard



120 6th Street • Fort Myers, FL 33907
T/239.489.0024 • Toll Free/800.306.0656 • F/800.306.0677
www.airscience.com

